

OTC HEARING AIDS USER MANUAL

Model No.: ZR-E101



Keep Health Always Around you!

Dear Customer

Thank you very much for choosing Our hearing aid.

Our hearing aids are high-quality devices designed specifically for people with hearing impairments. We are committed to the research, development, and production of high-quality hearing aids.

Please read the User Guide carefully before use and follow the correct methods for using and maintaining your hearing aid to extend its service life.

If you have any question, please be free to contact the local hearing aid service center or call our customer service center.

Basic Product Information

【Product Name】 In-ear hearing aid

【Product Specification and Model】 ZR-E101

【Product Dimensions】 25*20*15 millimeters


【Structural Components】 The in-ear hearing aid consists of a main unit, battery, and earpiece.

【Working Principle】 Sound signals are converted into continuously varying electrical signals via the microphone, then filtered and amplified by the circuitry. The receiver (headphone) restores these signals into amplified sound signals, which are then transmitted through the earpiece into the ear of the hearing-impaired patient.

【Indications】 After fitting, it is suitable for hearing compensation in patients with air-conduction hearing loss.

【Contraindications】 Acute otitis externa, otitis media, chronic suppurative otitis media (during the suppurative infection phase), acute suppurative otitis media, and individuals allergic to the materials used.

Main Structure

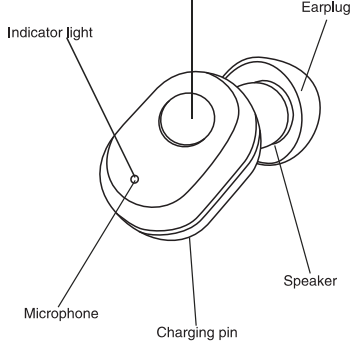

Product Specifications


Normal Working Conditions	
Ambient temperature	0℃~40℃
Relative humidity	≤85%
Atmospheric pressure	86kPa~106kPa
Working voltage	3.6V, 3.7V

Product Electroacoustic Performance	
Maximum OSPL90	<115dB-SPL
Output sound pressure level (SPL)	110dB~140dB
Frequency response range (Hz)	125~8000
Charge time with full charge (hours)	3.0~3.5
Equivalent input noise level	<20dB-SPL
Frequency response range (Hz)	125~8000
Frequency response range (dB)	40dB~100dB
Total harmonic distortion	<1.5%
Microphone battery	Low ohmic (silver zinc battery)
Capacity	15mAh, 3.6V/3.7V
Rated working current	<10mA

Installation Instructions


Choose earplugs of the appropriate size (there are several different sizes of earplugs in the box, choose the most suitable size for use to ensure greater comfort for your ears and better noise isolation).

Connect the earplugs to the hearing aids

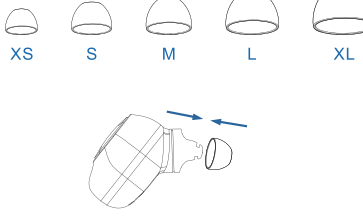


Operating Instructions

1

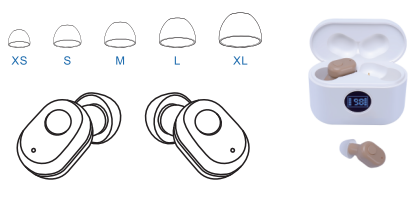


Remove the OTC hearing aids from the charging case. Select the earbuds that fit you from the five different sizes provided, connect the machine's connection port to the earbuds, and insert the earbuds into the inner ring.

2

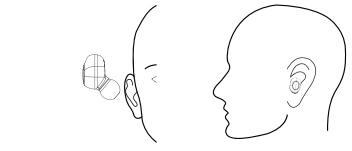


Slowly push the hearing aid with earplugs into the ear canal until it is close to the ear canal and feels comfortable.

3

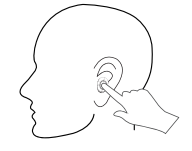

Press and hold for 3 seconds to activate the hearing aids (applicable to both left and right channels).

4

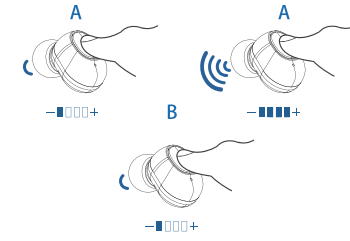

Touch the center of the casing to adjust the volume. There are 4 Volume levels to choose from:

A. Press the central button briefly and you will hear a single "beep," indicating volume level 1; press again and you will hear two "beeps," indicating level 2; the next level is three "beeps," indicating level 3; and four "beeps" indicate level 4 (maximum volume).

B. After hearing "beep beep beep beep" four times, touch the device again and the hearing aid will emit a single "beep" sound, restoring the volume to level 1.

When adjusting, do not cover the microphone hole with your hand, as this may cause feedback.

Hear The world Again

5

The above operations will be performed in sequence. After removing the hearing aid, please place it in the charging case to charge it promptly. The hearing aid will automatically shut down and begin charging. When the charging case is depleted of power, you must hold the button for 3 seconds to manually shut down the hearing aid; otherwise, it may cause feedback.



6

After placing the hearing aid in the charging case, the LED digital display will light up, and the orange indicator light on the hearing aid will also illuminate, indicating that the charging process is in progress. When the orange indicator light turns off, it indicates that the hearing aid is fully charged.

7



The LED display on the charging Case shows the charging status. When the charging case is charging via the connection cable, the numbers on LED display will light up and flash; when the number displays 100, it means that the charging case is fully charged.

Precautions, Warnings And Suggestive Instructions

- Below fitting hearing aids, you should go through professional inspections and hearing tests, and use them under the guidance of doctors or professional hearing aids.
- Moisture-proof, shock-proof and high-temperature-proof.
- Do not immerse the hearing aid in any liquid.
- Remove the hearing aid and put it in the charging case before swimming and wearing your hat and bath, after swimming or bathing, wait until the ear canal is dry before putting on the hearing aid.
- After the hearing aid is removed at night, it should be charged in a charging box, avoiding high light, high temperature, and high humidity.
- The outer surface of the hearing aid should be kept clean, and the eardrum in the earplug should be cleaned frequently.
- The functions should be used correctly, try to handle with care and avoid excessive force and excessive force.
- Hearing aids should be sent to the fitting center or dealer for maintenance and inspection regularly to ensure that the hearing aids are in good condition.
- This product should not be used for a long time at excessive volume to prevent damage to hearing. If you find discomfort, please stop using it.
- If the hearing aid is not used for a long time, the battery power may be discharged below the working voltage by the internal resistance of the battery. If the hearing aid cannot be turned on, please make sure that it is fully charged before trying to turn it on.
- If you hear a whistling sound in the ear, please check whether the size of the earplug (silicone/Melior) fits the user's ear canal and plug it tightly to prevent air leakage.
- It is not suitable to fall or collide, which may cause damage to internal components.
- It should not be opened without authorization.
- Hearing aids should avoid electromagnetic radiation, and do not use them in high-frequency thermal therapy or radiotherapy and strong magnetic fields.
- Never use it in a potentially explosive environment or in a dangerous state.
- Hearing aids cannot be treated as ordinary municipal garbage. Because hearing aids have built-in batteries, please use battery handling facilities to protect the environment.
- Do not let children play with this product and accessories, so as to avoid injury to children or possible damage to the machine.

Product Maintenance

- When the hearing aid is not in use, please charge it in a timely manner to ensure sufficient power for use.
- Clean the hearing aid casing with a clean, soft cloth.
- Clean and disinfect the earplugs with pure warm water.

Troubleshooting

Please follow the steps below to troubleshoot:

- Please confirm that the hearing aid can be turned on (the orange indicator light flashes 3 seconds and then goes out, indicating that it is turned on; the white version).
- Adjust the volume repeatedly to confirm whether the hearing aid produces sound.
- If there is no sound after the above two steps, ensure the device is fully charged before turning it on; then repeat steps 1 and 2.
- If there is still no sound after the above three steps, it can be determined that the hearing aid is not functioning properly or has malfunctioned. Please contact your local dealer or return the product to our company.

Transportation And Storage Methods

- Transportation Conditions: Handle with care to prevent drops; avoid direct sunlight, rain, or mechanical impact; keep dry.
- Storage Conditions: The product should be stored in a cool, dry place. The storage environment must meet the following conditions: (a) Ambient temperature: -10°C to 55°C; (b) Relative humidity: <85%; (c) Atmospheric pressure: 50 kPa to 155 kPa.
- Handle with care, avoid collisions or strong impacts.

Disposal Of Waste And Residue

This product contains a built-in battery and the main unit contains electronic components. Dispose of it in accordance with applicable national laws and regulations.

Package Contents

Hearing aids: 2

Portable mobile charging Case: 1 set

Earplugs: 2 each of extra small, small, medium, large, and extra large sizes

Cleaning brush: 1

User Manual: 1 copy

Name of manufacturer: SHENZHEN SHENHAI MEDICAL DEVICE CO., LTD.

Shenzhen, Guangdong, CHINA

Tel: +86-0755-89823193

Fax: +86-0755-81665233

Address: 8602 Building 1, No. 106 ZhongGuo Road, Xixiefu Community, Xixiefu Street, Bantou District 101, Shenzhen, Guangdong, 518000, CHINA